



PANalytical
get insight

Empyrean

Technical specifications



Enclosure



Specifications

Exterior dimensions	1400 (w) × 1162 (d) × 1947 (h) mm. The instrument is on wheels. It can pass a 989 mm door opening, the absolute minimum door opening is 867 mm.
Accessibility of interior area	1360 (w) × 1100 (h), allowing access to the experimental area by at least four people simultaneously (ideal for teaching)
Weight	1050 kg (2310 lbs.)
X-ray safety	Less than 1 microSievert per hour, measured at 10 cm distance, determined with a Mo source @ 60 kV, 50 mA
Safety loops	XSAFE: two independent safety loops, monitoring the conditions under which the shutter can be opened. XSAFE is set up in such a way, that operation of the instrument with open doors and open shutters is not possible, even for maintenance engineers. MSAFE: Best-in-class protection of users against unexpectedly moving parts by a double mechatronic safety interlock system
Cooling water supply	4 - 6 liters/min, pressure up to 8 bar, temperature 15 to 35 °C > dew point
Compressed air supply	House line, compressor or air bottle; 2-5 bar (0.2 - 0.5 MPa)
Power supply	Single phase 180 - 253 V; 47 - 63 Hz
Maximum power consumption (without controllers for optional equipment)	4.6 kVA

X-ray generation

Specifications	
X-ray generator	4 kW (max 60kV, max 100 mA)
Tube voltage	15 – 60 kV
Tube current	5 - 60 mA
Available anode materials	Cu, Co, Cr, Fe, Mo, Mn, Ag
Focus size	0.4 mm × 12 mm (LFF), 2 mm × 12 mm (BF)
Line/point focus rotation	Standard feature



PIXcel^{3D}



Specifications	
Detector size	256 × 256 pixels
Pixel size	55 µm by 55 µm
Point spread function	1 pixel
Count rate linearity	97% linear count rate per square mm: 13 million photons/(second mm ²)
Background noise	Less than 0.5 counts per second for the whole detector
Dynamical range	Up to 10 ¹⁰
Energy discrimination	Two-level discriminator for rejection of sample fluorescence, higher harmonics and white radiation (user adjustable)
Calibrations needed by user	None
Exchange of detection medium needed	None: solid-state detection technology

Goniometer

Specifications	
Configurations	Vertical goniometer, theta-theta and alpha-1 geometry
Measuring circle diameter (radius)	480 mm (r 240 mm) (can be changed for specific applications)
Goniometer angular range (w/out accessories)	360 degrees
Maximum usable range (depending on accessories)	$-111 < 2\theta < 168$ degrees
Angle positioning	DC motor drives with the next generation direct optical position sensing (DOPS2) using Heidenhain encoders and path tracking technology
Smallest addressable increment	0.0001 degree
2theta linearity over whole range	Equal or better than ± 0.01 degree
Maximum angular speed	15 deg/sec



Cradles

5-axes cradle (chi, phi, X, Y, Z)

Holds samples up to 80 mm diameter and 16 mm height, maximum weight 0.5 kg. With mapping range of 54 mm × 54 mm, or 2" wafers. The cradle has a chi tilt range of 96 degrees and phi rotation range of 720 degrees.

Options:


- dial gauge
- automatic height optimization
- microscope
- sample holder 80 mm (w or w/o levelling option)
- in-plane sample holder
- clamping device for PW18xx circular sample holders
- holder for 4×32 mm samples
- beam knife for line detector
- beam knife for reflectometry

3-axes cradle (chi, phi, Z)

Holds samples up to 140 mm diameter and 64 mm height, maximum weight 2 kg. Can handle non-ambient chambers like the DHS1100. The cradle has a chi tilt range of 96 degrees and phi rotation range of 720 degrees.

Options:

- sample holder 140 mm (w or w/o levelling option)
- transmission capabilities up to 40 degrees 2theta
- dial gauge
- automatic height optimization
- microscope
- clamping device for PW18xx circular sample holders
- beam knife for line detector
- beam knife for reflectometry

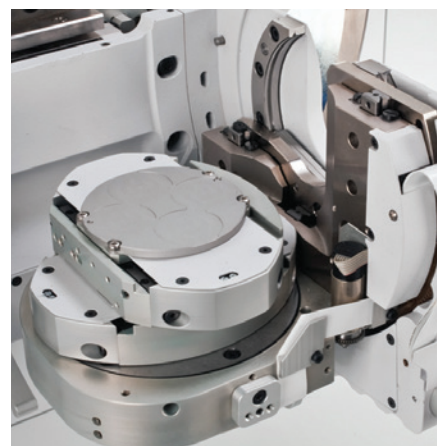

Empyrean

Cutting-edge technology in every aspect

In order to make the ultimate X-ray diffractometer for powders, thin films, nanomaterials and 3D objects, all essential components have been newly developed by our experienced R&D team.

- 1** Enclosure is on wheels for easy installation and relocation.




- 2** PANalytical's Empyrean Tubes offer robust line/point focus exchange.



- 3** The world's most accurate high-resolution goniometer

Minimum step size 0.0001° , 2θ linearity $\pm 0.01^\circ$

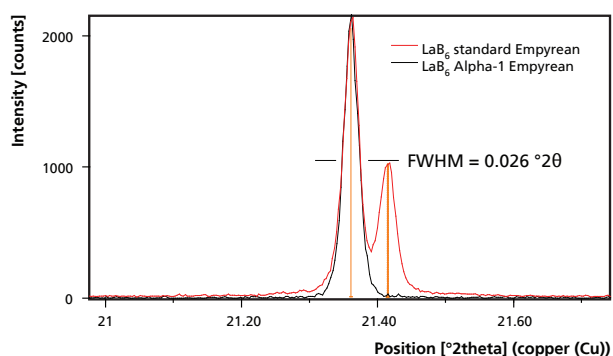

Empyrean sets the new bar for laboratory X-ray system resolution, with a FWHM of 0.026 degrees 2θ for the first reflection of NIST SRM660a LaB_6 .





- 4** The unique PIXcel^{3D} detector, co-developed with CERN and other leading European scientific institutes - more on next pages





5 New, unique PreFIX sample stages

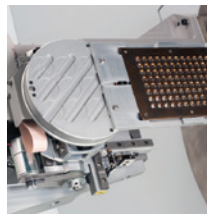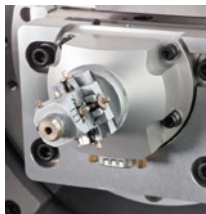

6 Universal PreFIX

The Empyrean system is equipped with universal PreFIX optics, stages and accessories, PANalytical's proven proprietary kinematic mounting concept for pre-aligned fast interchangeable X-ray modules.

7 Wide opening cabinet gives easy access to the experimental area.



8 Comprehensive software suite for data collection and analysis.



9 Meets all relevant worldwide regulations, for X-ray, electrical and mechanical safety.



10 19 inch rack mounts offering easy access to non-ambient controllers and integrated vacuum system

PreFIX optics

PreFIX (**pre**-aligned, fast interchangeable **X**-ray) is PANalytical's established method for robust, alignment-free component exchange. It allows for best performance, for every sample and application.

Incident beam optics, line focus

- Programmable divergence slit module
- Fixed divergence slit module
- Parabolic mirror (Cu, Co, others upon request)
- Focusing mirror (Cu, Co, Mo, Ag, others upon request)
- Hybrid monochromator 2xGe(220) for Cu, Co
- Hybrid monochromator 4xGe(220) for Cu, Co
- 4xGe(220) monochromator line focus asym. for Cu*
- 4xGe(220) monochromator line focus sym. for Cu*
- 4xGe(440) monochromator line focus sym. for Cu*
- Symmetric incident beam Johansson monochromator for Cu or Co


* Can also be used for point focus applications

Incident beam optics, point focus

- Monocapillary (available exit diameters from 50 μm to 1000 μm)
- Crossed slits collimator
- Double crossed slits collimator
- X-ray polycapillary lens
- 4xGe(220) monochromator line focus asym. for Cu
- 4xGe(220) monochromator line focus sym. for Cu
- 4xGe(440) monochromator line focus sym. for Cu


Diffracted beam optics

- Set of fixed receiving slits
- Set of fixed anti-scatter slits
- Programmable receiving slit
- Programmable anti-scatter slit with two independent motors, also for X'Celerator and PIXcel^{3D}
- Parallel plate collimators (0.09, 0.18, 0.27 degree)
- PreFIX 3-bounce sym. analyzer for Cu
- PreFIX 2-bounce asym. analyzer for Cu
- Rocking curve optics
- Diffracted beam parabolic mirror (Cu, Co)
- Diffracted beam monochromators for all detectors, including X'Celerator and PIXcel^{3D}
- 2nd diffracted beam path for automatic exchange between combinations of two diffracted beam optical modules
- Radius reduction interfaces 85 mm and 185 mm

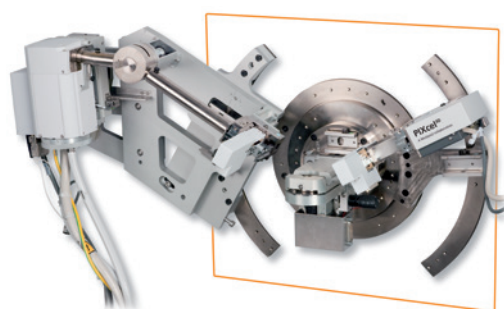

Detectors

	0D	1D	2D	3D
Xe proportional detector (optional monochromator)	✓			
Scintillation detector (optional monochromator)	✓			
X'Celerator (optional monochromator)	✓	✓		
PIXcel ^{3D} (optional monochromator)	✓	✓	✓	✓

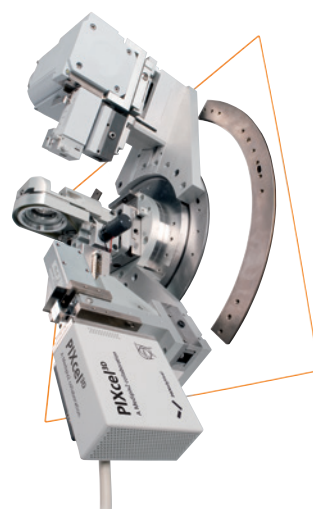

PreFIX sample stages

PreFIX sample stages can be easily removed from the goniometer and exchanged without realignment.

Non-spinning stages

- Programmable XYZ stage
- High-throughput stage
- Multi-purpose sample stage (MPSS)
- Pre-aligned capillary holder stage
- Flow cell for in situ crystallization research

Spinning stages

- Reflection transmission spinner stage
- Capillary spinner
- Microdiffraction spinner
- Computed tomography rotating stage

Stages with tilt capability

Empyrean's tilting stages are hot-swappable: they can be PreFIX exchanged and stored elsewhere without power interruption

- 3-axes cradle (chi, phi, Z)
- 5-axes cradle (chi, phi, X, Y, Z)

Sample stages for high temperature

All non-ambient chambers with Z-axis feature computer-controlled sample height optimization and automatic compensation for temperature-induced displacements of the sample

- HTK1200N rigid with Z-axis
- HTK1200N spinner with Z-axis
- HTK16N with Z-axis
- HTK2000N with Z-axis
- HTK2000N W with Z-axis
- Domed Hot Stage DHS1100

Sample stages for low and medium temperatures

All non-ambient chambers with Z-axis feature computer-controlled sample height optimization and automatic compensation for temperature-induced displacements of the sample

- TTK450 with Z-axis
- Oxford PheniX cryostat with Z-axis
- Cryo & Humidity Chamber (CHC Plus) with Z-axis

Transmission non-ambient stages

- VGI 2000M controlled temperature and humidity cell
- Cryostream Plus for cooling and heating of samples in glass capillaries



Reaction and high pressure stages

All non-ambient chambers with Z-axis feature computer-controlled sample height optimization and automatic compensation for temperature-induced displacements of the sample

- Reaction chamber XRK900 rigid with Z-axis
- Reaction chamber XRK900 spinner with Z-axis
- High Pressure Chamber HPC 900 with Z-axis

Sample changer for reflection transmission spinner stage

- 15 sample positions, upgradeable to 45 sample positions

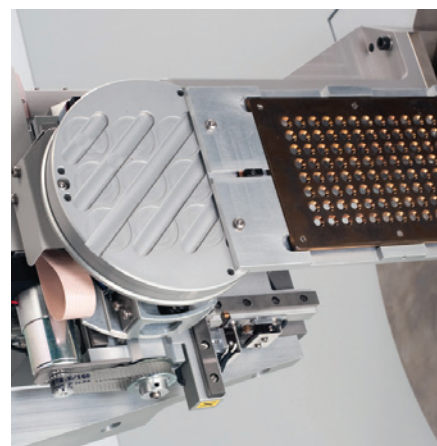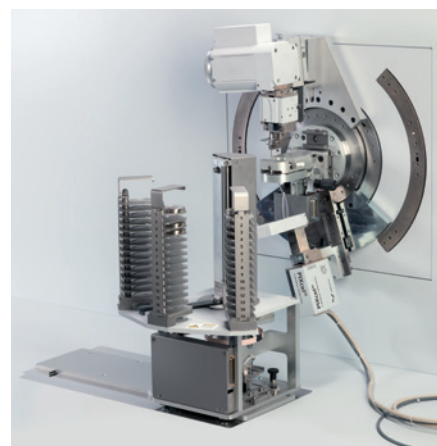

About PANalytical

PANalytical's mission is to enable people to get valuable insight into their materials and processes. Our customers can be found in virtually every industry segment, from building materials to pharmaceuticals and from metals and mining to nanomaterials. The combination of our software and instrumentation, based on X-ray diffraction (XRD), X-ray fluorescence (XRF), near-infrared (NIR) and optical emission (OES) spectroscopy and pulsed fast thermal neutron activation (PFTNA), provides our customers with highly reliable and robust elemental and structural information on their materials and is applied in scientific research and industrial process and quality control.

PANalytical employs over 1,000 people worldwide. The company's headquarters are in Almelo, the Netherlands. Fully equipped application laboratories are established in Japan, China, the US, Brazil, and the Netherlands. PANalytical's research activities are based in Almelo (NL) and on the campus of the University of Sussex in Brighton (UK). Supply and competence centers are located on two sites in the Netherlands: Almelo (X-ray instruments) and Eindhoven (X-ray tubes), in Nottingham, UK (XRF applications and standards), in Quebec, Canada (fusion sample preparation) and in Boulder CO, US (near-infrared instruments).

PANalytical is active in all but a few countries of the world. This worldwide sales and service network ensures unrivalled levels of customer support.

The company is certified in accordance with ISO 9001 and ISO 14001.

Visit www.panalytical.com for more information about our activities.

PANalytical is part of Spectris plc, the productivity-enhancing instrumentation and controls company.

Access to expertise

With the largest service network we are able to offer the most comprehensive support package possible.

Expertise:


- On-site training available
- XRF training courses
- Performance optimization
- Customizable expertise programs
- Assistance with multi-laboratory standardization

Care Agreements

Our customer support solutions have been developed with your business in mind. They are formulated as a family of four Care Agreements which can be tailored to your specific needs and provide fast, secure and reliable support.

- **ECONOMY:** indispensable coverage for self-sufficient operations
- **ADVANCED:** cost-effective support for routine usage
- **PREMIUM:** flexible package for high equipment usage
- **ELITE:** most comprehensive package for demanding environments

Global and near

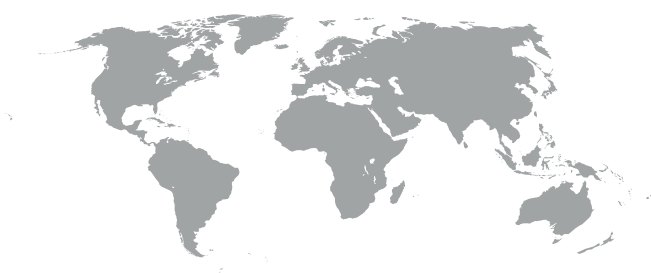

PANalytical B.V.
Lelyweg 1, 7602 EA Almelo
P.O. Box 13, 7600 AA Almelo
The Netherlands

T +31 546 534 444
F +31 546 534 598

info@panalytical.com
www.panalytical.com

Regional sales offices
Americas
T +1 508 647 1100
F +1 508 647 1115

Europe, Middle East, Africa
T +31 546 834 444
F +31 546 834 969

Asia Pacific
T +65 6741 2868
F +65 6741 2166